

EXTENSIONS AND ACCESSORIES

CORTEX SPIROERGOMETRY





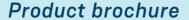

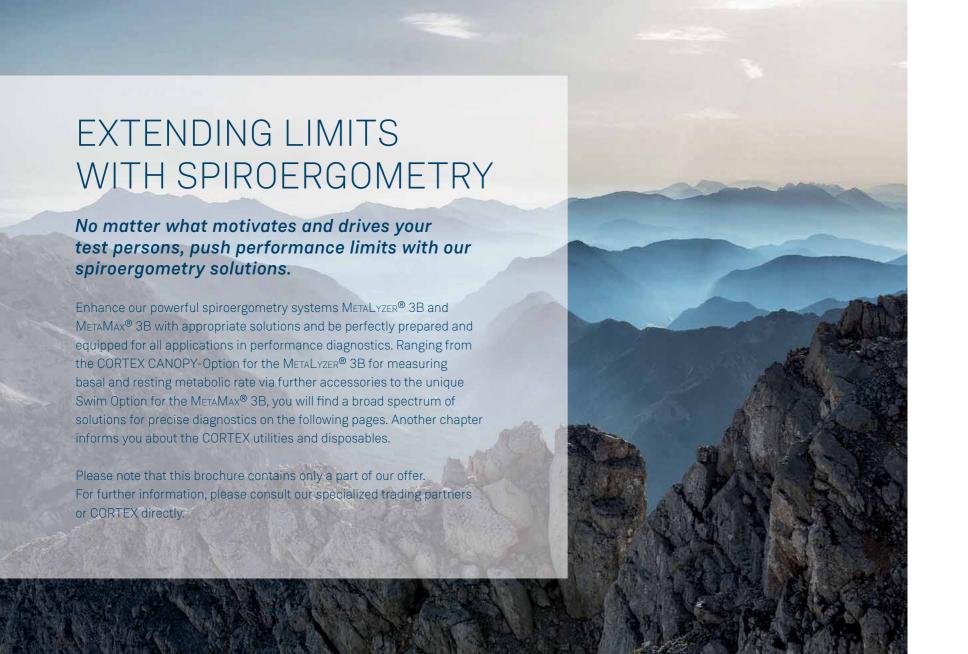

SYSTEM EXTENSIONS

Spirometry with the CORTEX
Spirometry kit

4-in-1 with the
Test Control Center

Swim Option for MetaMax® 3B

Measurement of resting metabolic rate with the CORTEX CANOPY-Option for MetaLyzer® 3B

from page 4



CORTEX RECOMMENDATIONS

Benefits through
Performance Packages

New in our assortment, the CORTEX disposable turbine

Ergometer
CORTEX Bike M

from page 12



SPIROERGOMETRY AND MORE

Cardiac output

Heart rate

ECG

Blood pressure

from page 18

SpO₂ measurement



MORE ACCESSORIES

Calibration kit

Sample Flow Dryer

Reusable turbine

Battery option

Equipment trolley

from page 24

SOFTWARE OPTIONS

Performance Assessment

Exercise Scheduler

Energy Metabolism

Spirometry

6-Minute-Walk-Test

Pre-operative
Risk Assessment

Network Licenses

HL7 Interface

from page 28



SPIROMETRY KIT



Spirometric examinations with an ergonomic mouthpiece

Expand your diagnosis with spirometric examinations using the practical spirometry kit and our hygienic disposable filters.

- no need to disinfect the volume transducer
- convenient handling
- compatible with MeтaLyzer® 3B and MeтaMax® 3B
- convenient, simple reordering system for disposable filters
- uncomplicated spirometric solution as a supplement to the spiroergometry

IAIN PRODUCTS AND ONSUMABLES	PART-NO.
and holder	010-00-100
SS Software Option Spirometry-kit incl. oftware license, hand holder and filter	606-01-001
SS Software Option Spirometry-kit incl. and holder and filter, without license	606-01-010
Safety-Pack II with in-moulded outhpiece, 1 carton à 100 pieces	500-00-029
ose clip	500-00-021
oftware license MSS Spirometry	721-09-028

TEST CONTROL CENTER (TCC)

Perform up to 4 spiroergometry tests simultaneously


The Test Control Center allows you to perform and monitor up to four spiroergometry tests at the same time.

Moreover, you will be completely free and flexible to choose which spiroergometry devices, ergometers and peripherals you want to combine with each other.

For example, you can use up to four Metalyzer® 3B,

Metal Max® 3B or any combination of the CORTEX devices and then combine these with treadmills or bicycle ergometers of your choice.

Each test can be started, monitored and stopped completely independently of the others, so that all test persons can pass their individual test. The software clearly presents all intuitively utilizable controls for the four tests.

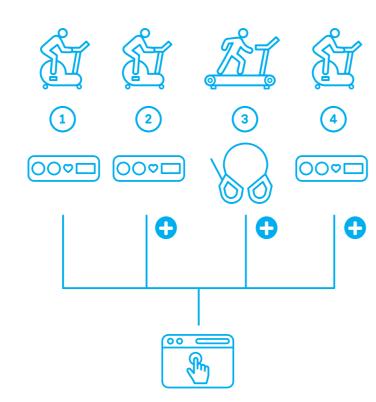

The Test Control Center is ideal for intensive training centers, sports centers and rehabilitation facilities, where a large number of athletes and patients are supposed to be looked after, or where many comparative tests are created within a narrow time window, e.g. during aptitude tests for young cadres.

Your benefits

- intuitive user interface in MetaSoft® Studio
- perform up to 4 tests simultaneously
- one operator can monitor up to four tests
- only 1 computer is necessary
- can be combined with MetaMax® 3B and MetaLyzer® 3B as required



TCC basic package

The basic package is your introduction to the Test Control Center (TCC). This includes the software license for MSS Test Control Center, a complete Metalyzer® 3B package with corresponding accessories and the application software MetaSoft® Studio, including connections, adapters and masks, as well as one BT Smart HR belt for measuring heart rate.

TCC extension

With the TCC extension package you can gradually purchase up to three additional CORTEX spiroergometry devices including the accessories required for the number of test persons.

MAIN PRODUCTS	PART-NO.
TCC basic package MetaLyzer® 3B (R3)	236-01-410
TCC extension MetaLyzer® 3B (R3)	236-01-420
МетаМах [®] 3B (R2) package	136-01-702

SWIM OPTION FOR METAMAX® 3B



The spiroergometry solution for swimming*

Benefit from the performance diagnostics through spiroergometry also in swimming. With the specially developed Swim Option for MetaMax® 3B, CORTEX enables exact performance physiological examinations even in water.

- suitable for a flow channel, pool and open water
- ergonomically designed special snorkels for freestyle and backstroke swimming with an integrated saliva suction device and splash-water protection
- integrated heart rate module with 5 kHz technology under water and Bluetooth® communication over water
- waterproof, shockproof transport case
- 12-hour operating time using battery option

THE SYSTEM CONTAINS ET AL.

Carrying case

Retrieval device, portable tube/cable holder

Battery option for MetaMax® 3B

Sample-Flow-Dryer, special sample line and dry granules

CORTEX-HR module for swimming, wireless charger

BT Smart HR belt

Freestyle and backstroke snorkels with splash-water protection

Volume transducer, disposable turbines, nose clips

MAIN PRODUCTS AND CONSUMABLES	PART-NO.
CORTEX Swim Option	632-01-010
Dry granules	010-00-023
Hydrophobic filter for special sample line	010-15-161
Special sample line (Tube in Tube), 3 m	632-01-080
Freestyle/backstroke snorkel	632-14-004
Backstroke snorkel	632-14-005
Disposable turbines (1 carton)	010-15-500
Sample-Flow-Dryer	132-01-081
CORTEX HR module swim	632-20-025
MetaMax® 3B	**

* Note: use only with the METAMAX® 3B (R2)



^{**} Note: detailed information can be found in our brochure MetaMax® 3B

CANOPY-OPTION FOR METALYZER® 3B



The extension for measuring the basal and resting metabolic rate

The CORTEX CANOPY-Option* for Metalyzer® 3B makes it possible to measure the basal metabolic rate and the resting metabolic rate without the test persons being affected by the mask.

- reliable measurement according to the gold standard
- particularly suitable for long measurements (without a mask)
- easy to clean and to disinfect
- fully automatic control of the hood ventilation
- high patient safety due to the automatic alarm function

THE SYSTEM CONTAINS ET AL.

Basic unit
Hood with veil
Tube
Bacterial- and viralfilter disposable (Pulmosafe)

MAIN PRODUCTS AND CONSUMABLES

CUNSUMABLES	PART-NR.
CORTEX CANOPY option (total package)	636-01-250
CANOPY-hood	636-01-001
CANOPY-veil for hood	636-01-015
CANOPY-tube	636-00-013
Pulmosafe bacterial- & viralfilter	500-00-019
MetaLyzer® 3B	**

DADT ND

* Note: use only with the MetaLyzer® 3B (R3)



 10

^{**} Note: detailed information can be found in our brochure METALYZER® 3B

CORTEXRECOMMENDATIONS

Do you already know our new, worldwide unique disposable turbine for spiroergometry? – The hygienic solution!

To make it easier for you to choose from our constantly growing product and service portfolio, we offer you the most popular combinations as pre-configured packages.

DISPOSABLE TURBINE



The first disposable turbine worldwide for spiroergometry

The singly packaged CORTEX disposable turbines in the practical dispenser box saves you the trouble of cleaning, disinfecting and calibrating a conventional turbine.

This will considerably reduce the amount of work you do between two measuring sessions, and all this takes place under the best possible hygienic conditions.

- no cleaning and disinfection necessary
- perfect hygiene
- simple handling easy to change
- no calibration necessary
- starting from 100 pieces in a convenient dispenser box
- suitable for both athletes and patients

TECHNICAL DATA

Peak flow	up to 14 l/s
Accuracy	±3% or 50 ml (whichever is the greate
nside diameter	30 mm

MAIN PRODUCTS AND

ONSUMABLES	PART-NO.
DRTEX disposable turbine carton à 100 pieces*	010-15-500
ask adapter for CORTEX flow sensor th disposable turbine	010-00-101
ectronic for volume transducer and sposable turbine MetaLyzer [®] 3B	236-01-021
ectronic for volume transducer and sposable turbine MetaMax® 3B	136-01-021

^{*} Note: bulk prices for 1000 or 5000 turbines

The Service Package for high investment and planning security

The Service Package contains a stock of consumables and accessories that are meant to enable one to work as independently as possible during the first years. The warranty extension of 3 years* which is also included in the package, gives you additional planning and cost security for your projects and studies. In addition, the operating costs are significantly reduced by the package.

* Note: this can further be extended as required, see warranty extension




Performance Package MetaLyzer® 3B

The Metalyzer® 3B Performance Package already includes the most commonly used hardware and software options in sports and preventative medicine. In addition, the package contains one calibration kit, one mask extension set and a larger stock of consumables.

Performance Package MetaMax® 3B

This package includes products contained in the Performance Package Metalyzer® 3B as well as other accessories specific to the MetaMax® 3B, such as a GPS receiver, back carrying system and a stationary kit for the optimal use of the MetaMax® 3B also in the laboratory with a 2-meter long sample line and a power supply unit.

ACKAGE	METALYZER® 3B	METAMAX® 3B
ervice Package	236-06-905	136-06-905
erformance Package	236-01-500	136-01-500

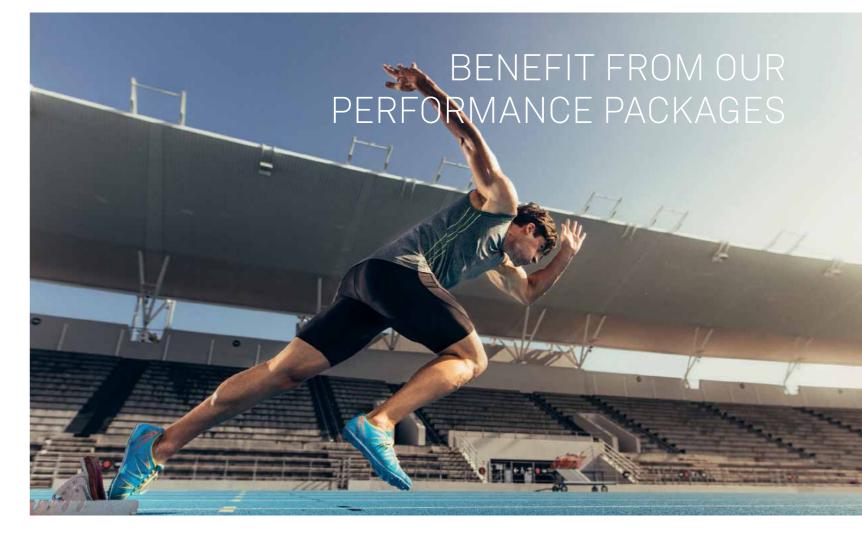

CORTEX WARRANTY EXTENSION

We provide a 2-year warranty on every CORTEX system.* With the CORTEX warranty extension, you can renew this time frame as you wish and without limitation.

CORTEX offers you attractive guarantee packages for the spiroergometry systems Metalyzer® 3B (R3) and Metalwar® 3B (R2), graded according to the scope and term. The protection covers the free repair and troubleshooting of hardware and software as well as regular, free software updates.

Take advantage of longer terms

Choose long terms as required and profit from an attractive price graduation.

Your benefits

- 100 % warranty directly from the manufacturer
- individual terms
- particularly attractive conditions for longer terms
- effective protection from unforeseeable repair costs
- free return of your device after repair



690-01-030
690-01-031
690-01-032

CORTEX BIKE M

The all-rounder among the ergometers

With its simple operating concept, low boarding height and electrically adjustable seat, the CORTEX Bike M is the ideal ergometer for daily use in clinics, hospitals and training facilities.

- open, low entry
- continuous, electrical seat height adjustment
- patient weight up to 200 kg
- high stability
- LCD display with remote control, ergometry and manual operation
- large workload range (6-1000 Watts)
- optional blood pressure measurement and/or SpO₂



MAIN PRODUCTS AND CONSUMABLES	PART-NO.
CORTEX Bike M	610-01-040
Automatic blood pressure measurement	610-05-051
oxygen saturation measurement	610-05-050
riathlon handlebar	610-00-011
Horizontal saddle adjustment	610-00-012
Quick release skewers without saddle	610-00-097

^{*} Note: in connection with the service package of 3 years

SPIROERGOMETRY AND MORE

The Metalyzer® 3B and the MetaMax® 3B provide a very detailed picture of the cardiopulmonary efficiency.

We also offer a variety of solutions for heart rate, electrocardiography and blood pressure diagnosis to provide a broader focus. A new addition is the non-compliance-based and non-invasive cardiac output measurement.

Here you can find a small selection of suitable devices that are easily compatible with our systems. In addition, our MetaSoft® Studio application software provides a large number of interfaces for measuring devices of other suppliers. Please contact us for a list of all deliverable or compatible manufacturers and models.

3

PHYSIOFLOW® ENDURO



Non-invasive measurement of cardiac output

The PhysioFlow [®] Enduro is a high-resolution impedance-cardiography device. Use this device to continuously measure cardiac output, stroke volume and ejection fraction non-invasively, non-compliance-based and without influencing the spiroergometry.

- non-invasive determination of the cardiac output
- objective findings independent of patient cooperation at rest and under load
- portable, lightweight system
- accurate assessment of the haemodynamic condition
- measurement on treadmill and ergometer

CONSUMABLES

NSUMABLES	PART-NO.
rsioFlow [®] unit	681-01-100
sor kits for PhysioFlow® Enduro*, 200 pieces**	681-05-001
ent cable for PhysioFlow® Enduro	681-05-010
etooth [®] antenna for PhysioFlow [®] Enduro	681-05-011
ch for PhysioFlow [®] Enduro	681-05-012

- * Note: Please use only original electrodes, otherwise we shall not be able to guarantee reliable results.
- ** Note: Sensor kits are available in different packaging units.
 Please inquire from us.

BT SMART HR BELT



Heart rate measurement

We recommend the use of a convenient chest belt for monitoring the heart rate. The heart rate is transmitted automatically to the MetaMax® 3B or MetaLyzer® 3B by means of a stable and energy efficient Bluetooth® transmission.

- Bluetooth® Smart-technology and encrypted 5 kHz-signal transmission
- comprehensive connectivity to physical fitness apps and diagnostics devices
- comfortable grip by means of an elastic belt, soft materials and secure locking system
- integrated connection to the spiroergometry software

MAIN PRODUCTS	PART-NO.
BT Smart HR belt	000-45-030

CUSTO GUARD



1- or 3-channel ECG

Using the 1- or 3-channel ECG custo guard, we offer mobile electrocardiography measurement in the practical sports belt.

- energy-saving wireless transmission to MetaMax® 3B or MetaLyzer® 3B
- simultaneous display of spiroergometry and ECG in MetaSoft® Studio
- light weight and reliable grip, even with extreme physical activity
- waterproof

MAIN PRODUCTS	PART-NO.
3-channel ECG set	840-01-020
1-channel ECG set	840-01-021

CUSTO CARDIO 400



Resting/Stress Test ECG

An improved pacemaker detection, the new custo true wave® filters, a sampling rate of up to 32 kHz and the high functionality transform the custo cardio 400 into a 12-channel ECG of the latest generation.

- a sampling rate of up to 32 kHz
- data transmission via USB and Bluetooth®
- defibrillation protection for dielectric strength of 5000 V
- suction capacity of 6 steps, 0-280 mbar
- electrode check with automatic indication of quality

MAIN PRODUCTS	PART-NO.
custo cardio 400	350-01-025
Software option MSS meta control 3000	721-09-013

Note: Further accessories like batteries, electrodes, suction unit etc., available on request.

CUSTO CARDIO 300



Mobile 12-channel ECG

The mobile 12-channel resting/stress test ECG has integrated LEDs which give information about the signal quality.

- Sampling rate up to 32 kHz for higher resolution and improved pacemaker detection
- start/stop button for ECG recording
- battery runtime of 10 hours

MAIN PRODUCTS	PART-NO.
custo cardio 300 BT accu - ERG*	350-01-076
custo cardio 300 BT battery - ERG*	350-01-077
Patient cable for cardio 300 (clip) long	840-05-016
Software option MSS meta control 3000	721-09-013

* Note: including short patient cable



NONIN® PULSE OXIMETER



Determine oxygen repletion

The Nonin® pulse oximeter WristOx2® Model 3150 delivers very accurate measurements for all sorts of patient groups and attitudes. It makes it possible to carry out cardiological supervision of out-patients, wireless remote monitoring and overnight feasible studies, both in the hospital and in a private setting.

- a well-thought-out, intuitive operating concept
- Bluetooth® 2.0 technology with a range of up to 100 m
- measuring technology for use in clinical studies, among others
- high skin compatibility according to the guidelines (RoHS)
- extremely robust technology for intensive, long-lasting use in clinics and private households

MAIN PRODUCT	PART-NO.
Nonin® Pulse Oximeter	671-00-011

SUNTECH® TANGO M2



Automatic blood pressure measurement

The Tango M2 conceived for stress tests on the treadmill and ergometer was specifically developed to filter out signal interferences caused by movements of the test persons.

- display of the systolic reading in less than 15 seconds,
 or of the whole blood pressure reading within 30 seconds
- automatic communication with MetaMax® 3B and MetaLyzer® 3B
- close measuring intervals selectable at time-critical moments or patient emergency situations

MAIN PRODUCTS	PART-NO.
Suntech® Tango M2	421-01-001
QRS-Trigger-Software custo med license	720-00-400
CORTEX-custo Trigger Cable Tango M2	420-01-011

BLOOD PRESSURE MEASUREMENT OPTION



In conjunction with the CORTEX Bike M

You can supplement the blood pressure option* using the ergometer CORTEX Bike M. The delivery includes a standard cuff. Other cuff sizes can be obtained on request.

MAIN PRODUCT	PART-NO.
Automatic blood pressure measurement	610-05-051

^{*} Note: only in conjunction with the CORTEX Bike M





The complete calibration solution

A complete set for the calibration of gas and volume flow sensors of all CORTEX spiroergometry devices, meant to guarantee optimal functionality of the systems.

- calibration gas syringe, 3 l
- calibration gas bottle
- gas-on-demand valve with adapter
- carrying case, calibration gas bottle holder

MAIN PRODUCTS	PART-NO.
Calibration Kit Professional with carrying case	010-01-793
Calibration gas bottle	010-00-005
Calibration syringe 3 I	010-00-083
Gas-on-demand-valve	010-01-898

SAMPLE FLOW DRYER





Reducing humidity in the breathing air

The accuracy of measurement is influenced by high humidity such as the one experienced in tropical or cold regions and in applications in swimming. To guarantee accurate measurement also under these conditions, the Sample Flow Dryer dries the inhaled air with the help of dry granules before it is analyzed by the sensors.

MAIN PRODUCTS	PART-NO.
Sample-Flow-Dryer	136-01-081
ry granulate	010-00-023
Sample line for Sample Flow Dryer, 60 cm	136-01-026
Special sample line or Swim Option, 3 m	632-01-080

REUSABLE VOLUME TURBINE





The reusable alternative to the disposable turbine

Use the volume turbine for stable and reliable measurement results in daily medical care and sports or also under extreme environmental test conditions, both indoors and outdoors.

- conceived for patients, amateur athletes and top athletes
- precise nano-production for stable and reliable test results
- suitable for all common disinfectants
- designed for MetaMax® 3B and MetaLyzer® 3B

TECHNICAL DATA

Peak flow	≤ 20 l/s
Accuracy	± 2 %
Inside diameter	28 mm

IAIN PRODUCTS	PART-NO.
ectronic assembly for MetaMax® 3B (R2)	136-01-007
ectronic assembly for MetaLyzer® 3B (R3)	236-01-013
eve	010-15-146

BATTERY OPTION FOR METALYZER® 3B

Be independent of the power grid

With the battery option, you are completely self-sufficient from the power supply during operation. Simply leave the device switched on if you would like to transport it from patient to patient without having to restart and recalibrate it every time. The battery option also prevents you from losing measured data in case of power failure, saving you the need to re-prepare the device and patient for testing.

MAIN PRODUCTS PART-NO. METALYZER® 3B with battery option 236-01-720 Retrofit set for battery option On request

CORTEX MOBILE PLUS

All devices in a practical trolley

The CORTEX mobile PLUS trolley offers sufficient space for all the testing and computer accessories. Apart from the spiroergometry device, the option enables the conveniently storing of one or two TFT screens, a computer and all other necessary accessories.

- large drawer for small components and consumables with kink protection for cables
- cable conduit integrated into the device column
- calibrating bottle and pump fixed under the drawer
- dual or single monitor mount

MAIN PRODUCTS AND CONSUMABLES	PART-NO.
Trolley CORTEX mobile PLUS with 19' TFT dual monitor bracket *	821-01-030
Trolley CORTEX mobile PLUS with 19' TFT single monitor bracket *	821-01-031
Printer shelf	821-00-030
Tray for MetaMax® 3B	821-01-015

^{*} Note: different country-specific versions available

METASOFT® STUDIO SOFTWARE OPTIONS WENCES

Like no other application software, MetaSoft® Studio combines the complexity of spiroergometric examinations with an extremely intuitive user interface.

With its versatile software options, daily test routines can be designed even more comfortably. If the focus is on the fast and reliable analysis of the medical data, the possibilities of using the Wasserman interpretation help, among others, or pre-operative risk assessment are useful options.







V'C02, HR

Oxygen Uptake



MSS PERFORMANCE ASSESSMENT



Effective training with the help of accurate training areas

The software option CORTEX MSS Performance Assessment identifies optimal training areas using established indicators (among others VT1, VT2, VO2 max, HRmax). On this basis, individualized heart rate training areas can be determined for different load phases and training goals.

- determination of individual heart rate training areas
- individual zone models can be generated
- accurate determination of VT1 and VT2 (AT and RC)
- manual adjustment of thresholds and training zones possible

MSS EXERCISE SCHEDULER



Enter concrete training guidelines

The CORTEX MSS Exercise Scheduler transfers the results of the MSS Performance Assessment into a graphic training plan which is organized like a time table. Based on the previously established performance level, defined framework conditions like training goal, load increase, frequency and type of sports activity, training plans are generated in the micro, meso and macro cycle.

- clear training plans in the micro, meso or macro cycle
- reference training plans for numerous sports types, prevention and rehabilitation measures
- motivation through clear training plans
- designing of personal training plans

MSS ENERGY METABOLISM FOR CPET



Optimize fat and carbohydrate combustion

MSS Energy Metabolism measures the individual energy metabolism, the combustion of fats and carbohydrates with increasing loads, and establishes the optimal heart rate zones for improved economization of the metabolism.

- determination of the increasing energy metabolic rate with increasing load
- determination of fat and carbohydrate combustion
- determination of individual heart rate zones for optimal energy burning and training effects
- motivation of patients through clear evaluations and reporting functions

MSS BASAL AND RESTING METABOLIC RATE



Determine the exact energy metabolic rate individually

As a reliable basis for a controlled calorie intake, e.g. in the context of a recommended nutritional plan, an accurate and individual determination of the resting or basal metabolic rate is imperative. This option measures the individual energy metabolic rate in a secure, intuitive and valid way.

- determination of the individual basal and resting metabolic rate with a mask or the CANOPY
- separate representation of energy, carbohydrate and fat consumption
- motivation of patients through clear evaluations and reporting functions

MSS SPIROMETRY



Perform a pulmonary function test

All CORTEX spiroergometry systems can also be used as spirometers. The software option MSS Spirometry offers VC, FVC and MVV tests as well as numerous other measurement, calculation and comparison options.

MSS 6-MINUTE-WALK-TEST

The 6-minute-walk-test is an approved method for the evaluation of COPD as well as for therapy control in rehabilitation. With our software option the advantages of spiroergometry are added to the system. Do not only measure the cardiopulmonary reaction of the test persons directly during the 6 minute walk test, but also compare it with the results of classical spiroergometry.

MSS PRE-OPERATIVE RISK ASSESSMENT

This license contains a module in which the assessor is guided through a partially automatic interpretation of results which is especially designed with pre-operative risk assessment in mind.

Make easier decisions on operative risks from a patient and get more certainty for planned interventions.

MSS NETWORK LICENSE BASIS



The MSS network license basis allows you to install the software on as many workplaces as required with one common database. Up to 3 users can work with the software simultaneously.

MSS NETWORK-EXTENSION




The MSS Network Extension allows you to increase the number of simultaneous accesses to the network version of MetaSoft® Studio.

MSS HL7 INTERFACE

The MSS HL7 Interface option provides you with extensive configuration options for connecting to a hospital information system. Please note that special effort can be calculated for setting up the interface.

SOFTWARE OPTION	PART-NO.
MSS Performance Assessment	721-09-020
MSS Exercise Scheduler	721-09-021
MSS Energy Metabolism for CPET	721-09-022
MSS Basal and Resting Metabolic Rate	721-09-023
MSS Spirometry	721-09-028
MSS Network License Basic	721-09-011
MSS Network License Extension	721-09-008
MSS 6-Minute-Walk-Test	721-09-035
MSS HL7 Interface	721-09-029
MSS Pre-operative Risk Assessment	721-09-025

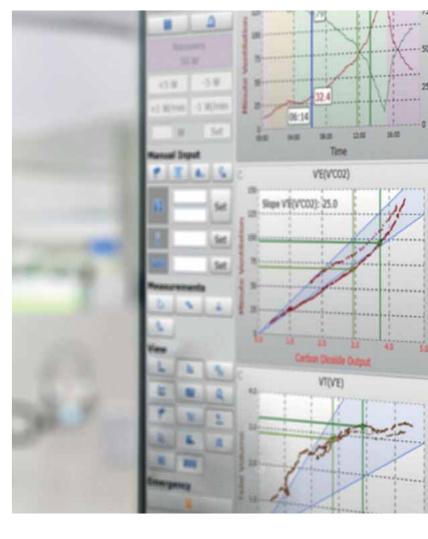

Volume Transducer	Hardware Options	Packages	Trolley & PC	Ergometer	HR/ECG	PhysioFlow/BD/SpO ₂	Software Options	Accessories & Spare Parts
CORTEX volume transducer disposable (p. 13) disposable turbines flow electronics CORTEX volume transducer reusable (p. 26) turbine flow volume electronic	battery option (p. 27) stationary kit (p.27)	performance package (p. 14) service package (p. 14) Swim Option for MetaMax® 3B (p. 8) Test Control Center (p. 6) for 2 CPET devices for 4 CPET devices warranty extension (p. 16)	CORTEX mobile Plus dual monitor (p. 27)	CORTEX Bike M (p. 17) - ③ seat height adjustment (motor) - ③ integrated BP - ③ integrated SpO2 - ③ adjustable pedal cranks - ③ triathlon handle bar - ③ adaptation for pediatric ergometry - ○ fall stop - ③ handrail - ○ TrackMaster - ③ Monark - ③ TrackMaster - ③ Micromed - ③ Marquette - ③ Quinton - ② Motek Force - ③ TechnoGym - ③	Polar H10 BT smart HF sensor (p. 20) MSS 1-lead guard (p. 20) MSS 3-lead guard (p. 20) custo cardio 400 (p. 21) MC3000 interface custo cardio 300 (p. 21) MC3000 interface ECG - electrodes custo med ECG - electrodes blue sensor 12-lead ECG interface Schiller GE Micromed AmedTec Dr. Schmitt Ergocid Spacelabs Norav Medset Labtech Welch Allyn	PhysioFlow (p. 19) PhysioFlow license PhysioFlow electrodes PhysioFlow cable bag batteries charger PhysioFlow electrodes (p. 19) for 50 tests for 200 tests Tango M2 for use with custo med (p. 23) QRS trigger software license Tango M2 with ECG module (p. 23) Nonin wrist 0x2 3150 (p. 22)	MSS Energy Metabolism for CPET (p. 30) MSS Performance Assessment (p. 29)	e face masks pediatric 1 pediatric 2 petite XS S M L e head cap adapter for face mask e sample line 0.6 m 2 m e turbine with sieve − reusable (p. 26) turbine − disposable (p. 13) 6 100 pieces 500 pieces 2 02 sensor Metalyzer® 3B 02 sensor Metalyzer® 3B 02 sensor Metalyzer® 3B bacterial filter spirometry PhysioFlow electrodes (p. 19) 6 for 50 tests 6 for 200 tests ECG electrodes custo med ECG electrodes blue sensor

The entire spectrum of spiroergometry

Configure your CORTEX spiroergometry system according to your individual requirements. You have a great variety of options and performance packages at your disposal.





The whole spectrum of spiroergometry

Visit our website www.cortex-medical.com and get to know the entire product range. Or ask for further brochures on other products. In addition to our mobile and portable spiroergometry solutions, we also offer customized solutions. Share your ideas with us! stamp:



IGZ Instruments AG Furtbachstrasse 17 8107 Buchs ZH

Tel. +41 44 456 33 33 igz.ch igz@igz.ch





Would you like to know more about CORTEX and our products?
Get in touch with us!
Tel.: +49 341 487 49-0 or on e-mail: sales@cortex-medical.com

<mark>ssign and specifications,</mark> misprints and/or errors excepted. Photos and diagrams are for illustra-<mark>are used in the docu</mark>ment. CORTEX Biophysik GmbH 2018.